

# Best Seller

## 2~8°C Pharmacy Refrigerator

High-quality Product Seven Sensors Temperature Uniformity Ensuring Safety of Storage



**YC-315EL**



Temp. Uniformity



Pulling Down to 5°C



Temp. Recovery After Door Opened for 1 Min

*\*Typical freezer data based on internal testing with freezer set point at 5°C and ambient temperature at 22°C.*

## Features



### Temperature Control System

Display temp.in 0.1°C increments



### Foaming Door



### Hydrocarbon Refrigerant

Environment-friendly  
Natural Gas



### CE Certificate



### Built-in Data Logger

Export up to 12 months of temperature data




### Alarms System

Multiple alarm functions

## Structure



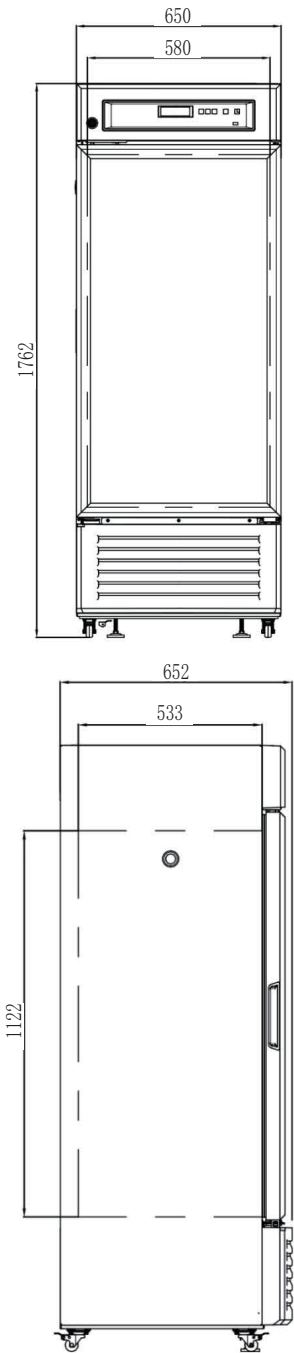

2~8°C Pharmacy Refrigerator		
Model	YC-315EL	
Capacity	315L/11.12Cu.Ft	
Internal Size (W*D*H) mm/Inch	580*533*1122/22.83*20.98*44.17	
External Size (W*D*H) mm/Inch	650*652*1762/25.59*25.67*69.37	
Package Size (W*D*H) mm/Inch	717*732*1815/28.23*28.82*71.46	
NW(Kgs/Lbs)	87/191.80	
GW(Kgs/Lbs)	99/218.26	
Performance		
Temperature Range	2~8°C	
Ambient Temperature	16~32°C	
Cooling Performance	5°C	
Climate Class	N	
Controller	Microprocessor	
Display	Digital Display	
Refrigeration		
Compressor	1pc	
Cooling Method	Air Cooling	
Defrost Mode	Automatical	
Refrigerant	R600a	
Insulation Thickness(mm)	D:55, L/R:35, B:52	
Construction		
External Material	PCM	
Inner Material	HIPS	
Shelves	5	
Access Test Port	1pc. Ø23mm	
Casters	4+(2 Leveling Feet)	
Data Logging/Export	Record every 10 mins/Max.12 months via USB port	
Interface	Remote Alarm Contact, RS485	
Backup Battery	Yes	
Alarms		
Temperature	High/Low Temperature, High Ambient Temperature, Condenser High Temperature Alarm	
Electrical	Power Failure, Low Battery, Sensor Error	
System	Data Communication Failure, Door Ajar, Data Logger Failure	
Electrical		
Power Supply(V/Hz)	110V/60Hz	220~240V/50Hz
Rated Current(A)	3.92A	1.35A

\*Zhongke Meiling reserves the right to change products and specifications without prior notice.